


# FEATURES & ASSEMBLY REFERENCE

LB2-PED-L2

**LARGE PEDIATRIC WRISTBAND**

After printing, push the wristband (where indicated) through the front and peel away from the back of the carrier sheet. Fold the laminate over to protect the image area and smooth all edges to create a tight seal.

Thread the end of the wristband through the slot. Size and secure by removing the wristband's adhesive tab.

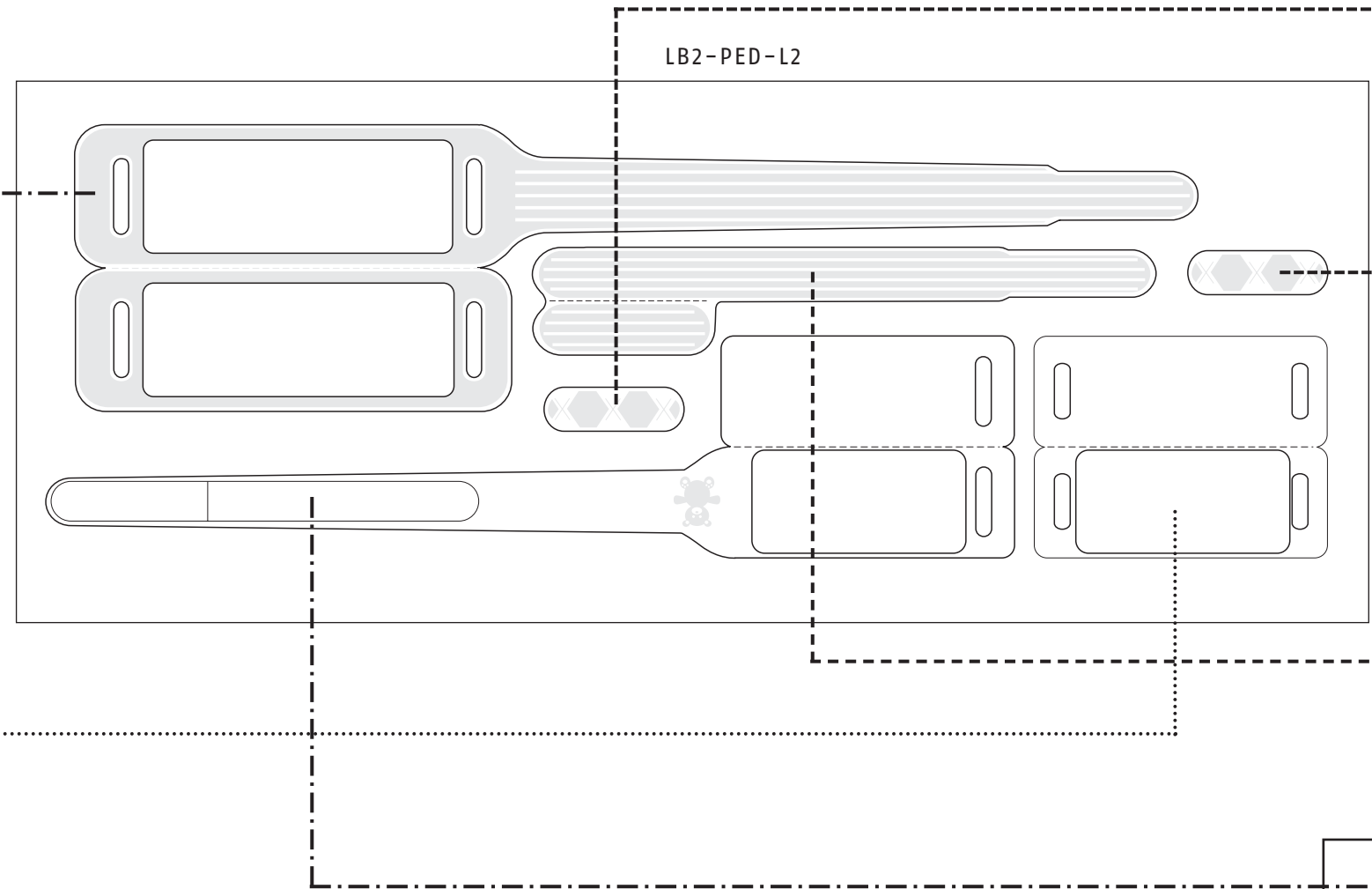

**SECURITY SEAL**

After the wristband has been applied to the patient, peel off the security seal located on the back of the form and pinch around the cinched end of the wristband to further secure the adhesive tab and indicate wristband tampering.

**WRISTBAND EXTENDER**

Simply push the extender (where indicated) through the front and peel away from the back of the carrier sheet. Line up the adhesive tab of the wristband with the extender, remove the wristband adhesive tab, and fold over the extender flap.

After wrapping around the patient's wrist, thread the end of the wristband through the slot. Size and secure by removing the extender's adhesive tab.

**PEDIATRIC ID TAG**

After printing, push the ID tag (where indicated) through the front and peel away from the back of the carrier sheet. Fold the laminate over to protect the image area and smooth all the edges to create a tight seal.

**COMFYCUFF® CONSTRUCTION**

Thread the end of the ComfyCuff® strap through the nearest slot, remove adhesive from strap, bond tag to adhesive, and thread end of strap through second slot. Wrap around patient's ankle or wrist, size, and secure by using the velcro.

For additional security and to indicate any wristband tampering, stick optional ComfyCuff® security seal across both ends of the velcro strap.

**SMALL PEDIATRIC WRISTBAND**

After printing, push the wristband (where indicated) through the front and peel away from the back of the carrier sheet. Fold the laminate over to protect the image area and smooth all edges to create a tight seal.

Remove only the adhesive tab section closest to the image area and begin wrapping around the patient's wrist. Thread the end of the wristband through the slot. Size and secure by removing the additional adhesive tab section. Trim any excess band.

LASERBAND RECOMMENDS LIMITING USE OF THE PEDIATRIC WRISTBAND TO THE PATIENT'S WRIST ONLY.



PRODUCT GUIDE


LB2-PED-L2

LB2-PED-L2

**PRINTING WRISTBANDS AND LABELS**

LaserBand® wristbands and labels have been designed to meet the rigid requirements of laser printers. All efforts have been made to assure that products leaving our factory have been manufactured as specified for sheet fed laser printers.

**LOADING INSTRUCTIONS:**

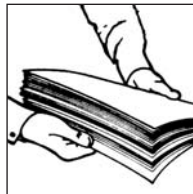


Remove the sheets from the protective plastic bag, fan the sheets on all four sides.




If a curl toward the wristband/label side is evident, bend the sheets before loading. Discard the top and bottom sheets, which may have a static charge.



Load sheets according to the instructions in your laser printer operator manual.

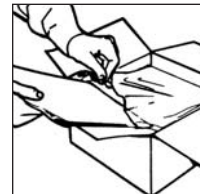

Upon completion of printing, return unused wristbands/label sheets to protective plastic bag.

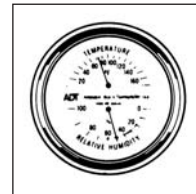

**STORAGE INSTRUCTIONS**  
Wristbands should be stored at 70°-80° F/50 ± 5% RH. Keep away from heat sources and do not expose to direct sunlight.


Wristbands and labels should be stored in the room in which they will be imprinted at least 48 hours prior to running for proper conditioning.

**FORMATTING AND LABEL DESIGN**

LaserBand® wristbands and labels have been designed to provide formatting flexibility. Care should be taken in the design stage to assure the desired copy will fit on the wristbands/labels.

The laser printers will hold a tolerance of ± 1/10". This dimension must be added to the dimensions of the desired imprint area to ensure a proper fit.

When bar codes or OCR codes are to be printed, an additional 1/10" must be added around the barcode to ensure proper scanning.

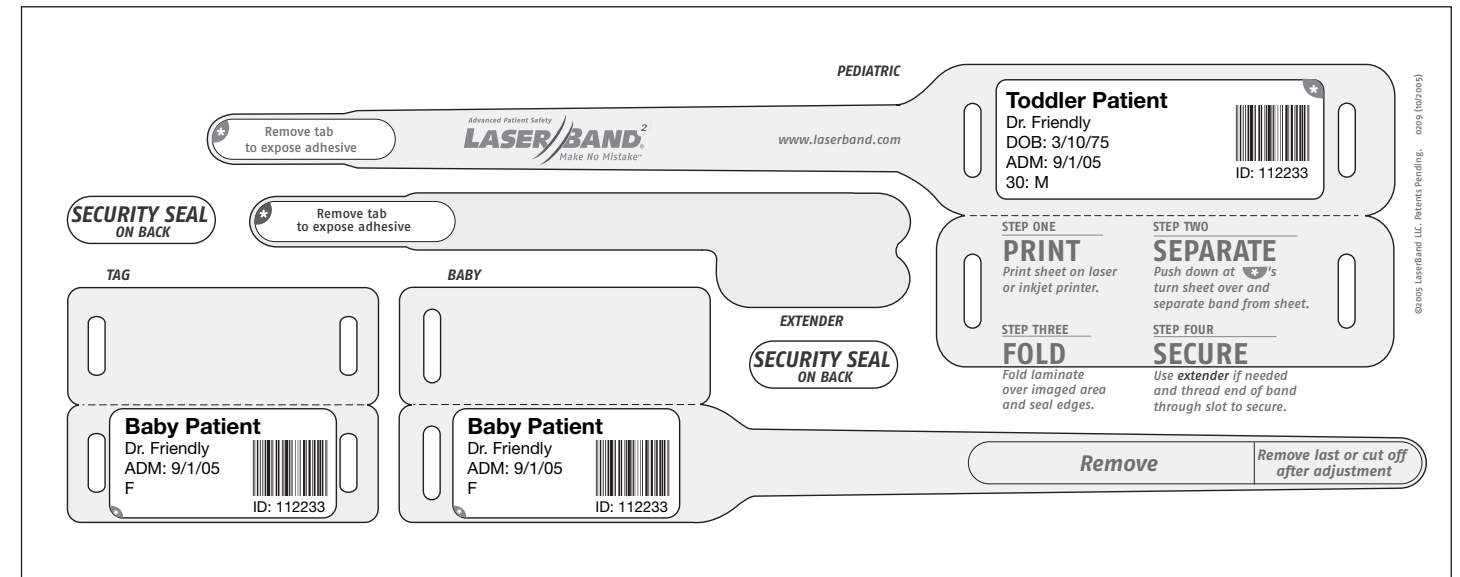

LOOK INSIDE FOR FEATURE AND ASSEMBLY REFERENCE.

